

DOCUMENTATION
for
PSYCHOPHYSIOLOGY PROTOCOL

in
MIDUS REFRESHER
BIOMARKER PROJECT
(P4)

University of Wisconsin ♦ Institute on Aging
(September 2017)

INTRODUCTION

This document provides an overview of the psychophysiology experimental protocol (laboratory challenge study) included in the MIDUS Refresher Biomarker Project (P4) protocol. It provides detailed information about the protocol and data processing procedures, as well as descriptions of variables created and basic guidance about their usage. Information is also included about construction and usage of administrative and computed variables.

Data users are also encouraged to review the Refresher Biomarker (P4) Readme Data File Notes. This document provides general information about naming conventions, as well as administrative and filter variables included in the data file. It also includes information about how we handled missing values and other issues that arose over the course of the study. For example, there are instances when variables were added or sections of an instrument were expanded for data entry purposes to accommodate additional information provided by the respondent.

This document will be periodically revised and updated as more information is gathered, and researchers continue to work with the MIDUS Biomarker data. If there are suggestions or comments, please contact midus_help@aging.wisc.edu.

TABLE OF CONTENTS

SECTION A: OVERVIEW OF DATA FILE AND COLLECTION PROTOCOLS.....	A-1
SECTION B: PSYCHOPHYSIOLOGY PROTOCOL FLOWSHEET and APPENDICES TO FLOWSHEET.....	B-1
SECTION C: DETAILED PSYCHOPHYSIOLOGY PROTOCOL DESCRIPTION.....	C-1
SECTION D: DATA QUALITY FILTER VARIABLES,,,.....	D-1
SECTION E: REFERENCES	E-1

SECTION A

OVERVIEW OF DATA FILE AND COLLECTION PROTOCOLS

OVERVIEW OF DATA FILE AND COLLECTION PROTOCOLS

The Biomarker Project (P4) Psychophysiology protocol is conducted in the morning on the second day of the CRU visit. The psychophysiology session is a standard, laboratory-based stress reactivity protocol and incorporates diverse measures from multiple sources as follows:

- Protocol Flowsheet (includes a Hand Usage Questionnaire)
- Saliva Samples for Cortisol assays
- Physiological measures of Stress Reactivity and Recovery
- Data Quality Filter variables

As described in the “Refresher Biomarker Project (P4) Readme Data File Notes”, the naming convention organizes variables according to the method used for data collection. We have followed this convention with respect to the psychophysiology data, thus analysts using saliva cortisol data in combination with flowsheet and/or cardiovascular reactivity data will need to pull variables from different sections of the data file as indicated below.

Protocol Flowsheet

The flowsheet variables appear in the data file immediately after the Actiwatch® data. Following the MIDUS Refresher naming convention, the variable names for the flowsheet data begin with the unique four characters set “RA4V”.

A copy of the flowsheet, with variable names added, can be found in Section B (below). Variable names generally appear to the right of the item they represent in brackets and bold capitals, however in some sections of the flowsheet they appear to the left of the item on the check line. See for example the second page of the flowsheet “Study start time [**RA4VST**]” or the fourth page [**RA4VS1T**] Collect Saliva Sample #1”.

The first page of the flowsheet, is a modification of the Edinburgh Handedness Inventory. This set of items was used to construct standardized indicators of laterality. Details about constructing laterality scores can be found in the “Documentation of Psychosocial Constructs and Composite Variables”.

The second page of the flowsheet, contains the session start time, questions about physical characteristics of the participant, and other factors that may influence experimental outcomes (e.g. consumption of caffeine, nicotine etc.), as well as a template providing an overview of the protocol order. The variables RA4VBEDH, RA4VBEDM, are new at the Refresher and represent the time that the participant got out of bed for the day. This time was added to the protocol to facilitate analysis of the saliva cortisol data.

New at Refresher: The time variables have been converted to a 24 hour clock with a restricted numeric format that allows leading zeros to be displayed. See ‘the Refresher Biomarker (P4) Readme Data File Notes’ for details about how time variables are handled in the date set.

This section of the data file also includes two flag variables. While cleaning the data we found a small number of cases for whom the session start time or the time that the participant got out of bed for the day was missing or later than expected due to a recording error. Thus, we imputed the times and then created the flag variables to identify the cases with an imputed value.

- RA4VSTF – flag for cases with an imputed Session Start time. The imputed time was computed using the site specific mode (50th percentile) of the lag time in minutes from the Session Start time to the 1st Saliva collection time where
 - Session time = 1st Saliva – Site specific lag in minutes
- RA4VWAKEF – flag for cases with an imputed Wake Time. The imputed time was computed using the average lag time in minutes from the blood collection time (recorded by CRU staff) to the Wake time where:
 - Wake time = Blood Collection time + Site specific lag in minutes.

The data file also includes computed two sets of computed lag variables indicating:

1. The number of hours from the time the participant last ate (RA4VATEH), last had caffeine (RA4VCAFH), or last had a cigarette (RA4VACIGH) to the Session Start time.
2. The number of minutes between Saliva sample collection times as well as the Wake up time, time out of Bed, and the Session start time.

The remainder of the flowsheet is the detailed protocol followed by staff during the psychophysiology session. Saliva samples were collected at designated points in the protocol, along with stress ratings. Staff recorded the saliva collection times and the stress ratings in designated spaces in the flowsheet for subsequent data entry.

Salivary Cortisol:

Saliva samples were collected for cortisol assay to provide a measure of neuroendocrine reactivity. Cortisol is a biomarker, thus these data appear in the data file with the other biomarker data immediately following the Physical Exam data. Consistent with the naming convention for the biomarker data, the variable names for the saliva cortisol values begin with the unique 4 character set "RA4B".

The saliva sample collection protocol was modified from MIDUS 2 to the Refresher in that 5, rather than 4 samples are obtained. This change was made to facilitate measurement of change in salivary cortisol in response to the stressors presented during the protocol. The samples were collected at five time points: 1) Pre-protocol (prior to attaching the ECG leads and other monitors); 2) Baseline (NEW): immediately before the Seated Baseline physiological recording begins; 3) Post-Cognitive Stress: after completing both Cognitive Stress Tasks and their corresponding Recovery periods; 4) Standing: immediately after the orthostatic challenge; 5) Recovery: 30 minutes after the orthostatic challenge. At the designated time respondents removed the cotton swab from the Salivette®, placed it in their mouth, rolled it around until saturated, and then put the swab back in the tube and replaced the cap. At the end of the session salivettes were stored in a -80°F freezer.

There is one other change from MIDUS 2, there is only one set of saliva cortisol variables, one for each of the 5 samples collected, and a final computed variable representing the average cortisol level across the individual samples. The saliva cortisol assay was run in duplicate for about half the MIDUS 2 Biomarker sample. We stopped running duplicates when it was determined that the cortisol assay provided high quality, reliable results and carried that practice forward into the Refresher.

For additional information about the saliva cortisol assay see the "Refresher Documentation for Blood, Urine and Saliva Data" and the "Refresher Biomarker (P4) Readme Data File Notes".

Psychophysiological Stress Reactivity and Recovery:

A detailed description of the psychophysiology protocol can be found in Section C. In addition to a detailed outline of the protocol it also provides descriptions of the measures and variables included in the data file, particularly the primary outcomes of interest, which are derived from the three physiological signals collected during the protocol: electrocardiogram (ECG), blood pressure (BP), and respiration (RSP).

The variable names for the stress reactivity data begin with the unique four characters set "RA4V".

This section of the data begins with two administrative variables:

- RA4ZPPHYS – categorical variable indicating whether the psychophysiology session was completed or not, and if not the point in the period in which the session was terminated.
- RA4ZPHYSD – categorical variable indicating whether the full set of physiological measures was obtained, and if not which measure is not available.

Data Quality Filter Variables:

Details about these variables can be found in Section D below. Two sets of filter variables are included in the data file. Although described at the end of this document, the first set of filter variables appear in the data file immediately after the above administrative variables (RA4ZPPHYS, RA4ZPHYSD). The second set appears between the HRV data and the Blood Pressure data.

SECTION B

PSYCHOPHYSIOLOGY PROTOCOL FLOWSHEET

And

Appendices to Flowsheet

FLOWSHEET

MIDUS Project 4 Psychophysiology Protocol

Date of Hand Usage Measure: ___/___/___
 mm dd yyyy

NOTE: Please administer Hand Usage measure the night before the psychophysiology session.

Hand Usage Questionnaire:

Please indicate your hand usage preferences in the following activities by circling the number in the appropriate column. If with any activity you use **both hands confidently**, mark the “Either hand or both hands” column. Some of the activities require both hands. In these cases, the part of the task or object for which hand preference is wanted is indicated in brackets. Try to answer all the questions, and only leave a blank if you have no experience at all with the object or activity.

	Strongly Left hand	Left hand	Either hand or Both hands	Right hand	Strongly Right hand
1. Writing [RA4VHAWR]	1	2	3	4	5
2. Drawing [RA4VHADW]	1	2	3	4	5
3. Throwing [RA4VHATH]	1	2	3	4	5
4. Scissors [RA4VHASC]	1	2	3	4	5
5. Toothbrush [RA4VHATO]	1	2	3	4	5
6. Knife (without fork) (e.g. cutting vegetables) [RA4VHAKN]	1	2	3	4	5
7. Spoon [RA4VHASP]	1	2	3	4	5
8. Broom [upper hand] [RA4VHABR]	1	2	3	4	5
9. Striking Match [match] [RA4VHASM]	1	2	3	4	5
10. Opening box [lid] [RA4VHAOB]	1	2	3	4	5
	Strongly Left foot	Left foot	Either foot	Right foot	Strongly Right foot
i. Which foot do you prefer to kick with? [RA4VHAKI]	1	2	3	4	5
	Strongly Left eye	Left eye	Either eye	Right eye	Strongly Right eye
ii. Which eye do you use when using only one? (e.g. using a camera) [RA4VHAOE]	1	2	3	4	5

Date of Psychophysiology session: ____/____/____
mm dd yyyy

Site ID:

UCLA 1
Wisconsin..... 2
Georgetown 3

Missing Data Codes: use throughout flowsheet

Don't Know = 7, 97, 997 etc.

Refused/Missing = 8, 98, 998 etc.

Inapplicable = 9, 99, 999 etc.

Study Start time [RA4VST] ____ : ____ a.m. p.m.

What time did you wake up for the day and not return to sleep?
[RA4VWAKE] ____ : ____ a.m. p.m.

What time did you get out of bed for the day? [RA4VBED] ____ : ____ a.m. p.m.

"What time did you last eat something?" [RA4VTEH]
Make sure participant had some calorie intake (juice/milk) before study. If not, provide snack to the participant ____ : ____ a.m. p.m.

"What time did you last drink a caffeinated beverage?" [RA4VCAFH] ____ : ____ a.m. p.m.

If a smoker, "What time was your last cigarette?" [RA4VCIGH] ____ : ____ a.m. p.m.

Participant Measurements: Height(cm) _____ Weight(kg) _____

Participant Age: ____ yrs DOB ____/____/____ (m/d/y)

Is Participant color blind? NO _____ YES [RA4VCLB] If Yes, blind to which colors? _____ [RA4VCLBT]

STAFF, next question is here in case of Finometer problems. May help understand a physiological reason for poor readings.

Does Participant have diagnosed Raynaud's Syndrome (circulatory disorder characterized by cold hands may affect Finometer BP data) or other diagnosed circulation problems? NO _____ YES [RA4VCIRC]

If YES, describe: _____

PROTOCOL SUMMARY TABLE

STUDY TEMPLATES: c:\Program Files\Ledona Solutions\63midusR_1.tpl OR 63midusR_2.tpl

GAcq files: C:\DATA\MIDUS\PHYSDATA **COM Ports:** GAcq-to-Actor= COM1 Acq PC-to-Finometer (finolink)=COM _____

Period #	Period Name	File Prefix	Period Description	P'pant Activity	Duration (sec)
1	C1	63?????C1	Waveform Calibration	Finometer setup	Wait
2	N1	63?????N1	NADA1	Signal check; Null height sensor; start BP	Wait
3	C2	63?????C2	Seated calibration	Arm BP; Resp spirotube	Wait
4	V2	63?????V2	Saliva Sample 2	Salivette	180
5	B1	63?????B1	Baseline 1	Sit quietly	660
6	M1	63?????M1	Stress Task 1: Math* (or Stroop)	MATHTurner task	360
7	R1	63?????R1	Recovery 1	Sit quietly	360
8	S1	63?????S1	Stress Task 2: Stroop* (or Math)	Stroop task	360
9	R2	63?????R2	Recovery 2	Sit quietly	360
10	V3	63?????V3	Saliva Sample 3	Salivette	180
11	N2	63?????N2	Standing Transition (NADA2)	Stand up; arm at heart height	Wait
12	C3	63?????C3	Standing Calibration	Arm BP; Resp spirotube	Wait
13	U1	63?????U1	Standing (Upright)	Stand still quietly	360
14	V4	63?????V4	Saliva Sample 4	Salivette	180
15	W1	Not saved	Regular end of protocol; or advance to T1 if needed	Wait 30 mins for saliva 5	Wait
16	T1	63?????T1	Protocol terminated early: high or low BP criteria met; monitor symptoms & BP	Sit or lie down quietly	360

* GAcq will automatically and randomly select order of presentation of Math and Stroop tasks.

Sampling Configuration: Summary

Signal	Amplifier	A/D Channel	Sampling rate (Hz)
ECG	ECG amplifier	0	500
BP	Finometer	1	500
Respiration (chest)	Inductotrace	4	20
Respiration (abdomen)	Inductotrace	5	20

Study Setup

Before participant arrives:

The day before psychophysiology session: instruct **P** not to drink any caffeinated beverages or be exposed to nicotine after midnight the night before the psychophysiology session.

Get MIDUS participant ID#; Get saliva cortisol supplies.

Turn on all the computers and physiological monitors. Check that all cable connections are in place.

Check X-keys Stick keypad on Stimulus PC. Run Macro Manager and confirm that the Stimulus PC “sees” the Stick Keys device and that the Associated Layout is: **C:\midus\Xkey\MIDUS_XKeyStick.xk8**

On Acq PC: Turn on the Finometer and start the Finolink software's Monitor function. Check that Finolink directory is set to **c:\data\midus\physdata**

Be sure data folder **C:\data\midus\physdata** on Acquisition PC is **empty** (all files from prior sessions should have been moved to individual folders named with Participant ID number).

Create new participant folder on Acquisition PC: In **C:\data\MIDUS\Archive[subID]** create a new folder with same Subject ID used for the Gacq output files: **63BCCCC** (B= site ID; C= MIDUS Subject ID)

After participant arrives:

When **P** arrives, interview him/her and fill in the demographic information on the first page of the protocol flowsheet.

Ask **P** about the Pacemaker item on the *P-phys Session Information Template*. Determine presence and operating schedule of any internal or external heart assistive device of any type. Common devices include: Left ventricular assist devices (LVAD; but can also be right-sided); cardiac pacemakers and defibrillators; vagal stimulators. Use notes space in the “If yes, what type” field to give as much detail as the participant can provide.

Ask **P** about history of fainting or issues with standing for long periods of time. Note any contraindications or safety concerns to running U1 period.

Inform **P** about the procedures and give instructions:

“In this part of the study, we are going to collect information on your heart rate, blood pressure and breathing while you rest quietly and during some challenging tasks. I am also going to place three electrocardiogram leads on you to measure your heart rate; two on your collarbones and one on your abdomen. Next, I am going to place stretch bands on your abdomen and your chest. These bands measure respiration. I will also place a blood pressure cuff on your upper arm and on your finger. During the lab session you will perform two challenging tasks and then you will move from the seated to the standing position. One task is a color-word matching task. The other is a simple arithmetic task. After you complete these two tasks, we'll ask you to get out of the chair and stand still for a few minutes. You'll be able to lean against the wall while standing.”

We monitor heart rate, blood pressure and breathing rate during all the activities I told you about. Whether or not we move to the next activity in the session depends on the readings at the previous stage. If any of your readings reaches a standard cutoff level specified for this study during any of the activities, I'll tell you that we will stop the session now. Then I will ask you to relax for a few minutes.

This is a standard procedure, so please do not worry if we stop the session. The data we collect will be very useful regardless of when the procedure is stopped and will help us understand better how stress affects the way the heart works in many different kinds of people.

Remind **P** that the session is about an hour and a half long, so if he/she needs to use the restroom, now is the best time.

Inform **P** that they may request to stop the session at any time if they feel extreme discomfort or cannot continue. Instruct **P** that if they wish to stop, please notify the technician.

_____ Ask **P** to change into hospital gown (if used at your site) or remove jackets or jewelry that interfere with being able to place ECG electrodes. Ask **P** to **turn off pager, cell phone, palm devices, and all other beeping devices**. Ask **P** to **remove wrist watch** and all jewelry from the arm from which you will collect blood pressure readings. Put **wristwatch out of view** so **P** will not be able to track clock time during the study.

Explain Stress Ratings:

During most of this test, you will be sitting quietly or performing the tasks. Periodically, I will ask you for a stress rating, which will be on the scale of 1-10 (with 1 being not stressed at all and 10 being extremely stressed). I will ask: 'may I have a stress rating please'? Then you will give me a number from 1-10 indicating your stress level at that given moment. Just give me the number. Don't elaborate.

V1 SALIVA 1 (Pre-protocol sample, STANDING) & Pre-Saliva1 stress rating.

Let's do the first stress rating now:

[RA4VPSSR] **"May I have a stress rating please."** _____



Ask **P** to stand up, then:

"The first thing I would like to do while you are standing is collect a saliva sample. Then, while you remain standing, I will place three electrodes on your upper body to collect heart rate data. I will also place respiration bands around your chest and abdomen to record your breathing."

Saliva Collection Instructions:

- a) Take cap off saliva collection tube (DO NOT Separate the two tubes), set aside
- b) **"Put the swab in your mouth and hold it there until I ask you to remove it – at least one minute -- until it is saturated. Do NOT chew or bite down on the cotton."**
- c) Put saturated swab back in tube; with the smaller tube inside the larger, outer tube; replace cap.

[RA4VS1T] **Collect saliva sample #1 Time _____ : _____ am pm (circle am/pm)**



Prepare Participant for physiological monitoring:

_____ Attach ECG electrodes and Inductotrace bands while **P** is standing; use alcohol preps to clean area.

_____ Attach Inductotrace cables to the Inductotrace unit.

1. White leads for Chest band, Black for Abdominal
2. Bands go outside of gown, but shouldn't restrict ECG cables

_____ Seat **P** in chair; encourage her/him to find a comfortable position.

"Please sit down in a comfortable position. Keep in mind that you will be sitting in this position for quite some time and therefore it is important that you are reasonably comfortable. I understand that it will be slightly uncomfortable with all this stuff on."

"Now I am going to place a blood pressure cuff on your non-dominant arm and on your middle finger. Please let me know if anything feels too tight. The finger cuff will pulse throughout the session. Every now and then I will relax it to give you a break, but your finger might become a little numb. If it feels too uncomfortable, please let me know."

Finometer Set-up:

- Wrap the arm cuff air hoses with the small loop in the center of the Velcro strap.
- Place arm cuff on upper arm at heart level on same hand as finger cuff; the label "artery" should be placed just above the inside of the elbow.
- Strap the frontend box on the forearm facing the ceiling, on **TOP** of the wrist.
- Wrap the Velcro strap around mid-forearm to secure the frontend box, cables and air hoses in position.
- Attach finger cuff cable to the frontend box; match red dots on the frontend receptacle and cable connector, and insert as far as it will go.

- Carefully insert air hose in other receptacle on frontend box; *plastic air hose is very fragile– handle with care!*
- Route the cable and air hose between two fingers to the frontend box.
- Wrap cuff on middle finger: point cable and tube toward wrist; center cuff between joints; cuff should cover both knuckles equally.
- Center LED (light emitting diode) and PC (photo cell), 2 dots on interior of cuff, symmetrically on sides of finger.
- Connect shoebox sensor to frontend box (a modular phone cord type outlet).
- Attach the pillbox height sensor at mid-armcuff or heart level.
- Attach the other square shaped sensor to the finger cuff.
- Place **P's** arm on the arm rest.

DEMONSTRATION AND INSTRUCTIONS

_____ Instruct **P** how respiration will be calibrated. Hand the respiration tube to **P**, then explain:

“Now we are going to practice calibrating the respiration monitor. You will be breathing in and out of this bag with your mouth six times. Make sure to breathe naturally and not forcefully, just enough to fill and empty the bag. First, inhale in a normal breath, then put your mouth tightly around the tube and exhale to fill the bag, then inhale to empty it. Keep your mouth around the tube and continue to breathe in and out for 6 full breaths. Be sure your lips are tight enough around the tube that no air escapes. Also, hold the tube at the top without touching the plastic bag. I will put this noseclip over your nose to make sure you breathe only through your mouth.”

Before having **P** practice, demonstrate the calibration procedure using your own respiration bag **and noseclip**. When finished, place noseclip on **P** and ask them to begin the calibration breaths. Count the six breaths out loud for **P**. At the end of each breath (one full exhale-inhale sequence), say **ONE, TWO**, etc.

_____ **Give explicit instructions for remainder of protocol, with the following points:**

1. Recording will be done mostly during quiet rest with a math task, a color-word matching task, and standing up.
2. The order of the math and color-word matching tasks is random and no one knows which test the computer will present first.
3. **The tests are designed to be difficult** and everyone makes mistakes during the session but **participants** should just keep going. They may feel a bit frustrated or upset at times

It is extremely important during this session that you refrain from moving as much as possible. Movement creates noisy signals from the electrodes, and if they are too noisy, we may have to restart the session. In addition, it is equally important that you do not speak during any of the tasks or resting periods except when I ask for a stress rating. Speaking out loud changes your respiration data and if it changes too often then we may have to start over. Also make sure not to cross your legs during the session since it affects heart rate.

Of course, if AT ANY TIME you feel sick, have pain, or there is anything that needs immediate attention, please speak up and let me know right away.

Do you have any questions?”

“Now I am going to move on to the tasks and briefly explain each one. Feel free to interrupt me and ask for clarifications when needed. “

_____ Place keypad in comfortable position relative to the dominant hand.

“Have you ever used a keypad like this one before? It’s similar to the keys on a computer or typewriter. Please familiarize yourself with the key pad. There are two sets of keys here that you will use, the colored keys on the right side, and the Yes-No keys near the middle. The other keys have no function – nothing happens if you press them.”

PRACTICE TASKS

_____ **Set up practice for the MathTurner task.** Give instructions **before** starting task on the Stimulus PC.

“We'd like you to perform a simple arithmetic task. The computer will show you a series of addition and subtraction problems. After the problem appears, you will see the word “equals”, then an answer to the problem will appear. Your task is to determine if the answer is correct or incorrect. If it is correct, you press “Yes” on the keypad. If it is incorrect, you press “No.” When the answer appears, you have only about one second to press Yes for “correct” or No for “incorrect,” then another problem will appear.

If you don't respond quickly enough, the computer will count your answer as wrong and will present another problem to you.

In this task, speed and accuracy are important. Concentrate on the problems and enter your answers as quickly as you can. Please do not speak at all during the task, and try to move as little as possible.”

“Do you have any questions?”

Answer whatever questions the subject may have. Emphasize the importance of not speaking during task.

“Let's try a practice session now.”

Click on “**Shortcut to MATHTurner practice**” on Stimulus computer. (Practice starts as soon as you run the Shortcut. There is no instructions screen as in the Stroop Practice.)

Practice MathTurner with keypad: Have the P do a practice session, but **NO MORE THAN 3** times, as needed; during the first trial observe to see where they are having trouble.

After practice: “**This task is designed to challenge you, so don't be discouraged if you make mistakes. Please concentrate and try as best as you can. Do you have any questions about the math task?**”

Answer any questions the P may have. Emphasize **importance of not speaking** during task.

Set up practice for the Stroop task. Click on “**Shortcut to Stroop practice**” on Stimulus computer.

“**The second task is a color-word matching task. The computer screen will present you with color names, for example, the word red or blue. These names will appear in different colors, that is, the word “blue” may appear in yellow letters (cue card). Your task is to press the key on the keypad which corresponds to the color of the letters. For example, if the word “blue” appears in yellow letters, you would press the key corresponding to “yellow”. There are four colors, as shown on the keypad: red, yellow, green and blue. Press the key that matches the color of the LETTERS in the word, not the color named by the word. During the task, the keyboard map of these colors will appear at the bottom of each screen.**

In this task, the computer will score your responses for speed and accuracy. If you don't respond quickly enough, it will score your response as incorrect and present a new problem.

Let's try a practice session.

Practice STROOP task: Start the practice session as shown on Stroop practice screen.

Do you have any questions about the color-word matching task?”

Answer whatever questions the P may have. Emphasize to **not speak and stay as still as possible** during task. Also, emphasize that in this task, speed and accuracy are important.

Standing task explanation: Explain that the final task simply is standing up for a few minutes. Explain that you (the researcher) will assist P in getting up out of the chair carefully. You will help P to lean against the wall, moving as little as possible. **Explain that they will be standing for about 10 minutes. Standing does not need to be practiced.**

SET UP PCs:

Stimulus PC: Set up ACTOR.

1. Double click on the ACTOR shortcut on the desktop.
2. Check that the Path setting in the Actor screen is exactly this: **c:\Progra~1\Ledona~1**
3. Click on "Listen for Requests." Then click on "Blank Screen" which causes the screen to go blank. The screen will stay blank until the stimulus programs are called or ESC key is pressed.

Acquisition PC: GAcq

1. Click on Shortcut to GACQ Template Toggler on desktop:
2. Enter SITE ID (UCLA = 1; Wisconsin = 2; Georgetown = 3)
3. Enter the SUBJECT ID. **GAcq will begin collecting data in the C1 period as soon as you enter Subject ID and press OK.**

Which GAcq template used? See top of GAcq window (circle one):

1

2



*****BEGIN DATA ACQUISITION*****

C1 WAIT CALIBRATION 1: FINOMETER SET UP

On FINOMETER screen, press the **[Calibrate waveform]** tab card button until you see the square wave graph appear as shown in the figure.



Activate the [Finometer-Research] instrument using the **"red"** configuration; this will automatically terminate the square wave calibration and start the data recording screen.

Enter participant data in the [Describe Subject] tab card. Use metric converter chart to put height in meters and weight in kilograms. **Remember to press [Describe Subject] button again or the data will not be saved!!**

ADVANCE TO NEXT PERIOD WHEN READY

N1 WAIT NADA 1

CHECK ECG and Respiration signals on Gacq screen.

If the ECG and/or respiration signals are poor, then halt data acquisition and move the electrodes and/or bands around. Then restart NADA 1 period.

[RA4VLC] Optimal ECG Lead Configuration (use & **circle** whichever one gives best signal):

1. Upper Right and Upper Left 1
2. Upper Right and Lower Left (just above ground) 2
3. Other (specify) 3

Move the upper electrodes closer to each other on a horizontal line (toward center of chest) and slightly lower if you get a bad signal with the standard arrangement.

FINOMETER START UP:

Prepare to start Finometer measurement. Warn the P: **"I'm starting the blood pressure monitor. You'll feel the finger cuff pulsing. Are you ready?"**

Null Height Sensor: Hold the shoebox and pillbox sensors at heart level and press the **"mark"** key to null height sensor, then place sensors in proper locations on finger and arm cuffs.

Start a measurement by pressing on the **"start/stop"** key.



Check Blood Pressure signal quality: You should see a display of BP signals on the Finometer screen AND on the Gacq screen.

Make sure all the signals are clean before proceeding.


1. If ECG and/or respiration signal is poor, move electrodes and/or bands around.
2. If blood pressure signal is poor, press "start/stop" key on Finometer and move the cuffs around.

ADVANCE TO NEXT PERIOD WHEN READY

C2	WAIT	CALIBRATION 2 (SEATED)
		Inform the P: Now I'm going to calibrate the blood pressure readings. The armcuff will now inflate and stay inflated for about 2 minutes. It will be slightly uncomfortable. Also ask them to remain still and quiet during the calibration.
		Calibrate brachial arm pressure:
		a. Press on the [Physiocal & RTF-cal] button. Look at the left side "Physiocal" column to be sure it is set "on" (by default it should be ON. If not, change setting).
		b. Use the right arrow to activate the "RTF-cal section" and then select the "step" inflation type by pressing the [Physiocal & RTF-cal] button.
		Calibrate Respitrace: hand the P the respiration tube and place the nose clip on him/her. Instruct him/her to "wait for the initial count before blowing in and out of the tube naturally" . Return to Acquisition PC and begin counting.
		<i>As soon as you begin the first count:</i>
		1. PRESS AND RELEASE F1 : WHEN SUBJECT BEGINS FILLING THE BAG WITH THE 1ST BREATH.
		2. PRESS AND RELEASE F1 : AT END OF LAST BREATH.
		Right before advancing to V2 period: TURN OFF PHYSIOCAL: Press [Physiocal/RTF] tab card button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.
		NOTE: When you press the [Physiocal/RTF] tab card button once, this may toggle the Physiocal between ON/OFF setting. You may not need to use the arrow keys to move to the Physiocal section.
		After calibration, you cannot move the respiration bands from their position until the next calibration.
		ADVANCE TO NEXT PERIOD WHEN READY

V2	180	SALIVA 2 (Baseline Saliva Sample, SEATED)
		"Put the swab in your mouth and hold it there until I ask you to remove it – at least one minute -- until it is saturated. Do NOT chew or bite down on the cotton."
[RA4VS2T]	160	Collect saliva sample #2. Time ____ : ____ am pm (circle am/pm) 
	30	TURN ON PHYSIOCAL: Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to ON.
	15	TURN OFF PHYSIOCAL: Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.
	10	"Please sit quietly and try to relax. Remember not to speak unless you need help."
B1	660	<u>BASELINE 1</u> Participant sitting quietly.
	30	TURN ON PHYSIOCAL: Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to ON.
[RA4VSRB1]	25	"May I have a stress rating please." BASELINE _____ 
	15	TURN OFF PHYSIOCAL: Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.
	10	"Please Sit quietly and remember not to speak. The first task is about to begin."

M1 or S1 360 **STRESS TASK #1:**

[RA4VCS1] TASK ORDER: which task is it? Stroop or MATHTurner (circle one). 

_____ 30 **TURN ON PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to ON.

[RA4VSRCS1] 25 “May I have a stress rating please.” TASK #1 _____ 

_____ 15 **TURN OFF PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.

_____ 10 **“Please sit quietly and try to relax.” (Wait til task finished to say this !)**

R1 360 **RECOVERY 1**

_____ 30 **TURN ON PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to ON.

[RA4VSRR1] 25 “May I have a stress rating please.” RECOVERY #1 _____ 

_____ 15 **TURN OFF PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.

_____ 10 **“Please sit quietly and remember not to speak. The next task is about to begin.”**

M1 or S1 360 **STRESS TASK #2:**

[RA4VCS2] TASK ORDER: which task is it? Stroop or MATHTurner (circle one). 

_____ 30 **TURN ON PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to ON.

[RA4VSRCS2] 25 “May I have a stress rating please.” TASK # 2 _____ 

_____ 15 **TURN OFF PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.

_____ 10 **“Please sit quietly and try to relax.”**

R2 360 **RECOVERY 2**

_____ 30 **TURN ON PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to ON.

[RA4VSRR2] 25 “May I have a stress rating please.” RECOVERY #2 _____ 

_____ 15 **TURN OFF PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.

V3 180 **SALIVA 3 (Post-stress and recovery, seated): Saliva Sample**

“Put the swab in your mouth and hold it there until I ask you to remove it – at least one minute -- until it is saturated. Do NOT chew or bite down on the cotton.”

[RA4VS3T] 120 COLLECT SALIVA SAMPLE #3 Time ____ : ____ am pm (circle am/pm) 

_____ 30 **TURN ON PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to ON.

_____ 15 **TURN OFF PHYSIOCAL:** Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.

N2 WAIT **NADA 2: STANDING TRANSITION**

“Now we will transition to the standing portion of the test.”

Help the **P** get into the standing position, making sure that the connections to all the cables are maintained. **Have PP rest arm with Finometer on flat, stable surface such as hospital table or countertop. Arm rest should be located on non-dominate side . Pillows may be used to prop arm higher if needed as long as arm is stable. Move chair behind PP.**

Set up an arm rest at Heart Height for **P** to rest the arm with Finometer cuffs attached. Adjust their arm so the HITE reading on Finometer is as close to 0 as possible **(height difference between pillbox and shoebox sensors).**

Instruct PP to notify technician if they feel faint, dizzy, or have extreme discomfort. Inform PP a chair is located behind him/her and to sit down if discomfort or light headiness occurs.

Check All Signals:

1. You should see a display of BP signals on the Finometer screen AND on the Gacq screen.
2. Check quality of ECG, BP and Respiration signals on the Gacq screen.
3. If the ECG signal and/or respiration signal is poor then move the electrodes and/or bands around.
4. If BP signal is poor, press “start/stop” key to stop Finometer recording. Move the cuffs around.

ADVANCE TO NEXT PERIOD WHEN READY

LOW BLOOD PRESSURE CRITERIA CHECK:

IF BP falls below 80/60 OR participant feels faint, dizzy, weak, nauseous, OR heart rate suddenly slows to below 60 beats per minute (displayed on Finometer device), and remains at that level for 5 seconds etc. *during NADA 2 or anytime after P is standing*, help participant to sit down in chair. If needed, have participant lie down or bend and place head between legs from seated position. Proceed to **TERMINATION period in GAcq and Termination page of the flowsheet.**

C3 WAIT **CALIBRATION 3 (STANDING)**

Inform the **P**: **Now I’m going to calibrate the blood pressure readings again like I did when you were seated. The armcuff will now inflate and stay inflated for about 2 minutes. It will be slightly uncomfortable.** Also ask them to remain still and quiet during the calibration.

After at least 100 seconds of standing, calibrate the arm blood pressure:

Calibrate brachial arm pressure:

- a. Turn **Physiocal ON** (leave on until instructed to turn off after resp cal below).
- b. Use the right arrow to activate the “RTF-cal section” and then select the “**step**” inflation type by pressing the [Physiocal & RTF-cal] button.

Calibrate Respitrace: hand the **P** the respiration tube and place the nose clip on him/her. Instruct him/her to “**wait for the initial count before blowing in and out of the tube naturally.**” Return to Acquisition PC and begin counting.

As soon as you begin the first count

1. **PRESS AND RELEASE F1:** WHEN SUBJECT BEGINS FILLING THE BAG WITH THE 1ST BREATH.
2. **PRESS AND RELEASE F1:** AT END OF LAST BREATH.

Place the P’pant’s hand with the finger cuff at heart height. Use adjustable height table, an arm sling or other method that is comfortable. ***Make sure the arm is at heart level.***


Check All Signals, and adjust anything if needed.

Leave Physiocal ON for at least 15 sec after P’pant is stabilized and all signals look good.

Right before advancing to U1 period: TURN OFF PHYSIOCAL. Press [Physiocal/RTF] tab card button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.

Once you complete calibration, you cannot move the respiration bands.

ADVANCE TO NEXT PERIOD WHEN READY

U1	360	<u>STANDING TASK</u>	
_____		[If BP falls below criterion level, or P feels faint, dizzy, weak, nauseous, etc., help P to sit down in chair. If needed, have P lie down or bend and place head between legs from seated position. Proceed to TERMINATION period below.]	
_____	30	TURN ON PHYSIOCAL: Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to ON.	
[RA4VSRUV1]	25	“May I have a stress rating please.” STANDING _____	
_____	15	TURN OFF PHYSIOCAL: Press [Physiocal/RTF] button. Use arrow keys to move to left section labeled Physiocal and change setting to OFF.	

V4 **180** **Saliva Sample 4 - In STANDING position**

“Put the swab in your mouth and hold it there until I ask you to remove it – at least one minute -- until it is saturated. Do NOT chew or bite down on the cotton.”

[RA4VS4T] COLLECT SALIVA SAMPLE #4 Time _____ : _____ am=1 pm=2



_____ Set a timer for 30 minutes as soon as Saliva 4 collected.

_____ Leave monitors and Gacq running through entire V4 saliva period with Physiocal OFF.

[RA4VTERMV4] If session was terminated early, after which period of the protocol was Saliva sample #4 collected?



Circle One:

C1 N1 C2 V2 B1 M1 R1 S1 R2 V3 N2 C3 U1 W1 T1

W1 **WAIT** **WAIT / PAUSE BETWEEN STANDING AND TERMINATION PERIODS**

After Standing task ends, Gacq advances to this period and COUNTS UP FROM ZERO indefinitely until you Halt Acquisition. This period is used to end a normal session and prevent Gacq from advancing automatically to the Termination period. The Termination period should not run unless you manually advance to it when needed. See below.

HALT Acquisition for a Regular Session: If the session ran normally (no early Termination), end the GAcq session at any time during this WAIT period.

1. Click the **Halt Acquisition** button in GAcq.
2. **Stop Finometer data collection:** Press on the “start/stop” key (but leave finometer power on for data transfer later),
3. **Power off the ECG machine.**

Saliva Sample 5 + **30 MINUTES** AFTER SAMPLE #4 – In SEATED position.

“Put the swab in your mouth and hold it there until I ask you to remove it – at least one minute -- until it is saturated. Do NOT chew or bite down on the cotton.”

[RA4VS5T] COLLECT SALIVA SAMPLE #5 Time _____ : _____ am=1 pm=2



NOTES AND COMMENTS ABOUT THE SESSION: (you can enter these in GAcq Session Notes in the Information Template doc. Label as Session Notes)

END OF REGULAR SESSION

GO TO NEXT SECTIONS OF FLOWSHEET FOR:

1. Termination period (if needed)
2. Termination Debriefing (if applicable)
3. Disconnecting participant from equipment
4. Transfer Data

T1 TERMINATION: END SESSION DUE TO BP CRITERION OR OTHER PROBLEMS

If the BP goes **above** (usu. in Baseline or a Task period) or **below** (usu. during Standing period) the criterion level, complete this section of Flowsheet. **Termination is a type of recovery period. The HPV technician must do several things in quick succession and often in parallel:**

_____ **Press F1 key when you first notice BP in criterion range, then press it AGAIN when full termination criteria are met (i.e. when you make decision to terminate session & just before you move Gacq to the Termination Period).**

Above-criterion BP level: SBP _____ DBP _____ GAcq period: _____ GAcq Counter: _____ 

_____ **“May I have a stress rating please.” RATE STRESS WHEN TERMINATION CRITERION REACHED.** 

1. Encourage P to sit quietly and relax. Goal is to see if BP returns to a range in the study criteria within 300 to 240 seconds on GAcq counter (1 to 2 mins) and to the P’s near-baseline BP levels by end of the Termination period.
2. Determine if Participant is experiencing symptoms associated with high BP (or any other discomfort).
3. Monitor BP on Finometer. Use chart below to record BP every 30 sec or so.

_____ If BP criterion is reached during a **task**, inform the P: **“Thank you. You may stop responding to the task now. I am going to stop the task on your screen”**. If criterion met during Baseline, proceed below.

ADVANCE GAcq to Termination period: Click on Termination period in Gacq. This should halt whatever task is running on the Stimulus PC, if applicable.

_____ Inform the P: **Thank you. You met the requirements for this part of the session. We are now starting the last phase of this session.**

_____ **Finometer BP at beginning of Termination:** _____ 

_____ **Symptom report:** Ask P how s/he feels, if s/he has any discomfort. If YES, probe for symptoms such as: dizzy, feel flushed, nauseous, any pain of any kind. To get symptom report, strike a balance – try not to alarm the P by prompting him/her too specifically about types of symptoms, but we need to know if the high BP reading is also accompanied by such symptoms.

Symptoms: _____

Inform the P: Please sit quietly now and try to relax. Keep your eyes open and do not let yourself doze. You can better judge P’s functional status with eyes open and awake. If they doze off, wake them gently and be sure they feel OK. **If s/he can’t seem to keep from dozing, ask a nurse or doctor to see P right away.**

Record BP reading from Finometer every 30 seconds, as Gacq counts down. Don’t worry if you don’t get every reading or don’t get them all exactly on time.

Counter: 330 BP: _____	Counter: 210 BP: _____	Counter: 90 BP: _____
Counter: 300 BP: _____	Counter: 180 BP: _____	Counter: 60 BP: _____
Counter: 270 BP: _____	Counter: 150 BP: _____	Counter: 30 BP: _____
Counter: 240 BP: _____	Counter: 120 BP: _____	Counter: 0 BP: _____

Continue monitoring BP until end of Termination in GAcq. If BP does not recover to near-baseline levels or P reports symptoms that do not remit by end of period, (dizzy, feel flushed, nauseous) have doctor/nurse see P.

_____ **Did you call a doctor or nurse?** No _____ Yes _____

If yes, why, and what was the outcome? _____

_____ 0 ***If BP returns to a range in the study criteria within 300 to 240 seconds*** on GAcq counter (1 to 2 mins) and to the P's near-baseline BP levels by end of the Termination period, proceed with the next steps:

Ask for a stress rating.

“May I have a stress rating please.”

RATE STRESS AT END OF TERMINATION PERIOD. _____



Inform P: This part of the GCRC visit is complete. Thank you very much. You met the criteria for this visit. Now I'm going to ask you to do 2 more saliva samples.

_____ **HALT Acquisition for a Terminated Session:**

1. Click the **Halt Acquisition** button in GAcq.
2. **Stop Finometer data collection:** Press on the “start/stop” key (but leave finometer power on for data transfer later),
3. **Power off the ECG machine.**

_____ **After Termination Period, collect Saliva Samples 4 and 5, if participant is feeling well enough. Use saliva sample section shown above after Standing Period section of flowsheet.**

END OF TERMINATION PERIOD

TERMINATION DEBRIEFING

If session was terminated early, at this point give P more specific information about why. If the BP went too high or too low (during standing), explain what happened and tell them their BP value which pushed them over the criterion. Explain that this does not necessarily mean that they have a medical problem, but that we recommend that all participants who reach that level during the session be advised to see their own physician to have their BP and heart checked.

Notes about Termination Period: (you can enter these in GAcq Session Notes in the Information Template doc. Label as Termination Notes)

DISCONNECT PARTICIPANT

_____ **Power off** the Inductotrace and ECG.

_____ **Halt Finometer data collection.** Press the “start/stop” key on the Finometer, then press hard on both double arrow keys simultaneously to stop measurement and return to the start display. Do NOT yet power off the Finometer.

_____ Remove Finometer height sensors, finger and arm cuff.

_____ Disconnect finger cuff electrical and pneumatic connections.

_____ Remove frontend from wrist.

_____ Remove all ECG leads, tapes and bands from the participant.

_____ Disconnect Inductotrace OUTPUT cables from unit (**if you don't, the battery will run down!!**)

_____ Hit ESC key on Stimulus computer and tell Actor to Stop Listening.

TRANSFER DATA

Psychophysiology Session Info File:

- a. Complete the session information form. Save the PDF file: 63BCCCCC_info.PDF
- b. Session information sheet should include any restarted periods and reasons for restarts. Also include reasons for signal quality issues, i.e. "respondent was restless during M1" Note any equipment or software malfunctions. Include steps taken to resolve signal quality issues.
- c. Save the file to: C:\data\MIDUS\subID]

Transfer Finometer blood pressure files from Finolink to Acq PC.

If the Finolink Monitor software was running during the session, the finolink file(s) should already be on the Acq PC in c:\data\midus\physdata. If you had to restart measurement on the finometer during the session, there will be more than one file.

If the Finolink Monitor software was NOT running during the session, move the file(s) now from Finometer to Acq PC (the order of these actions is important!) :

1. Start Finolink with shortcut on Acq PC's desktop. Choose **Download** option.
2. Select Configure menu, then Directories. Be sure both directories are set to: **C:\data\MIDUS\physdata**
3. Click **Connect** button on right. You should see a list of files stored on Finometer. Click on Date and/or Time headings to sort files so you see the most recent file at top. Click current session's file(s) to select it (there may be more than one file!!)
4. Look at **Local Files** in left panel. It should display **C:\data\MIDUS\physdata**. If not, click button and browse to choose this folder.
5. Click on the **< (left arrow)** button in middle of screen to send file from Finometer list (right) to Acq PC (left).

Power Off Finometer: Once Finolink file(s) is on the Acq PC, on the Finometer monitor, exit to the start display by pressing hard on both double arrow keys simultaneously. Switch OFF the Finometer with the switch at the rear.

At this point, ALL Gacq and Finometer files from this session should now be in **C:\data\MIDUS\physdata**.

Stimulus PC Data Cleanup:

Copy task output files from **Stimulus PC** to archive folder on **Acquisition PC**.

63BccccM1.MathTurner & 63BccccS1.STROOP → **c:\data\archive\SubID for this participant]**

(That is, each P'pant has a folder on Acq PC labeled with their Subject ID. Their Stroop, Math and all physio files will be stored in this folder.)

63BCCCCC (B= site ID; C= MIDUS Subject ID)

Move stress task output files to **archive folders** on the **Stimulus PC** in these folders:

Stroop (.stroop):* **c:\data\archive\Stroopperfdata**
MATHTurner (.MathTurner):* **c:\data\archive\MATHperfdata**

(That is, all Stroop files for all participants archived in one folder, all Math files in other folder)

Now, there is a copy of the task output files both in the participant's folder in Acq PC archive AND in their respective Task Archive subfolders on Stim PC.

Acquisition PC Data Cleanup:

_____ **Rename Finometer files:** Append Subject ID plus an "underscore" character to ***beginning*** of Finometer filename. In Windows Explorer, select file in list, then Right Click, select Rename, then type new name.
Example: If SubID is 63190001, change **2jl28p08.271** to **63190001_2jl28p08.271**.

_____ **MOVE** all data files from **C:\data\MIDUS\physdata** to **C:\data\Archive\SubID**.
The physdata folder ***should now be empty***.

_____ Now, **ALL** the participant's files -- Stroop, Math, Gacq and Finometer -- from this session should be in **C:\data\MIDUS\physdata**. **Total files for a normal session:**

5 files per period (0,1,4,5,EVT) X 14 periods; 1 finometer file; 1 Session Info PDF file; 2 task files; 1 EVT file for W1 period (data files not saved but evt file is saved).

TOTAL = 75 files

If extra Finometer files were generated, you will have more than 75 files. If session terminated early, less than 75 files.

_____ Backup *all* data from this session, all of **C:\data\Archive SubID** to a USB Flash Drive or other medium preferred by your site.

_____ Turn off the computers and the monitors.

File Transfer to Columbia (CUMC) & Flowsheet to UW:

_____ Send all files from the **c:\data\archive\SubID for this participant** folder, the same ones you backed up to the USB flash drive.

_____ Copy the data files from this session onto the UW FTP server under your site's folder

Be sure there is a ***subfolder named with the Subject ID*** and all files are inside that folder.

_____ Send original Flowsheet to Madison site for archiving and data entry. The MIDUS Madison team will enter all flowsheet data into data bases.

_____ Keep a photocopy of all flowsheets on file at your site as a backup

APPENDICES TO FLOWSHEET

APPENDIX A

MIDUS Biomarker Project Revised Psychophysiology Blood Pressure Termination Criteria December 13, 2005

The Finometer can produce brief artifactual values for one or two SBP-DBP cycles. These aberrant values, if they are really artifacts, should begin to correct to more normal levels within a few heartbeats/BP cycles. Thus, the criteria below include a 5-second duration at various points to allow staff to determine with greater confidence that they are “real” physiological values and the session should be terminated.

Termination Criteria for HIGH Blood Pressure

Baseline Exclusion Criteria:

Subjects will be excluded from participation if their baseline blood pressure is 180/100 (either systolic OR diastolic criterion is met) or greater. This will be based on 3 manual BP readings taken at least 5 minutes apart, with subject seated.

Session Termination Criteria:

Persistent “low” level increase - The session will be terminated if blood pressure rises above 200/110 (either systolic OR diastolic criterion is met) and persists at that level for 1 minute, or the respondent complains of chest pain, vision changes, and/or headache.

Immediate Termination - The session will be terminated **immediately** if blood systolic blood pressure rises above 210 and does not begin to fall within 5 seconds. That is, if ONLY systolic BP reaches 210, you terminate. If diastolic BP reaches greater than 110 but systolic stays below 200, the criteria specified at #1 above should be applied.

Termination Criteria for Low Blood Pressure

Immediate Termination Criteria:

The session will be terminated at any point during the protocol if either systolic or diastolic BP falls 20 mmHg compared to either average baseline levels or the average level in the preceding 5 minutes, AND remains at that lowered level for 5 seconds without **beginning** to rise again, while the subject’s monitored finger is held at heart level.

Sitting to Standing Transition Termination Criteria:

The session will be terminated during the Sitting to Standing transition if the above criteria are met or:

The BP falls below 80/60 (both systolic and diastolic fall below this criterion), or

The participant appears distressed or complains of feeling lightheaded or faint or clammy and the Finometer blood pressure is falling steadily, or

There is a sudden slowing of the heart rate to below 60 beats per minute (displayed on Finometer device), and remains at that level for 5 seconds

APPENDIX B

SALIVA COLLECTION PROTOCOL & EARLY TERMINATION OF PSYCHOPHYSIOLOGY PROTOCOL

Introduction

The numbering of the saliva samples is linked to when they are collected during the psychophysiology protocol. When the psychophysiology protocol is terminated early it is important that the salivette numbers and saliva collection times be recorded appropriately, so that any samples collected can be linked to corresponding samples from R's who completed the protocol as written.

There are three types of situations in which the psychophysiology protocol is not completed. These are described below along with the procedure to follow for collecting saliva samples and recording times.

1. The session is not run due to an equipment failure or, at the last minute the R becomes unable to participate in the protocol (i.e. becomes ill) and so the Baseline period of the protocol is never begun.
 - a. If saliva sample #1 has been collected, record the time and send to the Biocore as usual.
 - b. If saliva sample # 1 has not been collected, do NOT collect it.
 - c. Do NOT collect saliva samples 2, 3, 4, or 5.
2. R is terminated early due to high (or low) BP
 - a. During baseline period before beginning the cognitive challenges. *Ideally saliva samples 1, 2, & 5 are collected.*
 - i. If possible collect a sample 30 minutes after termination. Use tube #5, which corresponds to the final sample, collected 30 minutes after the protocol has ended. Record the time in the location for sample 5 collection.
 - ii. If a sample cannot be collected write "Early Termination – Not Collected" at the location for recording the sample 5 collection time
 - iii. Do NOT collect Sample #3 or Sample #4, write "Early Termination – Not Collected" at the location for recording the sample #3 and sample #4 collection times.
 - b. During one of the cognitive challenges. *Ideally saliva samples 1, 2, 3, & 5 are collected.*
 - i. If a sample can be collected use tube #3 which corresponds to the sample collected at the end of the Saliva 3 period. Record the time of sample collection at that location as well. (Samples 1, 2, 3, & 5 are collected).
 - ii. If a sample cannot be collected, write "Early Termination – Not Collected" at the location for recording the sample #3 collection time.
 - iii. Do NOT collect Sample # 4, write "Early Termination – Not Collected" at the location for recording the sample # 4 collection time.
 - iv. Collect Sample # 5 30 minutes after the 3rd sample, record the time in the location for the sample 5 collection time.
 - c. During the standing challenge. *Ideally all 5 saliva samples are collected.*
 - i. If a sample can be collected use tube #4 which corresponds to the sample collected at the end of standing period. Record the time of sample collection at that location as well.
 - ii. If a sample cannot be collected, write "Early Termination – Not Collected" at the location for recording the sample # 4 collection time.
 - iii. Collect sample # 5 as usual.
3. R is in a wheelchair or unable to stand for some reason and is therefore unable to do the Standing Challenge. *Ideally all 5 saliva samples are collected.*
 - a. Collect sample # 1 and sample # 2 & #3 as usual.
 - b. Terminate the protocol after Saliva 3 (during NADA 2).
 - c. Collect sample # 4 about 15 minutes after sample 3, record the time in the space for sample # 4 collection time.
 - d. Collect sample #5 30 minutes after the 4rd sample, record the location for the sample #5 collection time.

SECTION C

DETAILED PSYCHOPHYSIOLOGY PROTOCOL DESCRIPTION

Introduction

The MIDUS Biomarker Project (P4) psychophysiology session is a standard, laboratory- based stress reactivity protocol. The data were collected at UCLA, Georgetown, and the University of Wisconsin and processed at the Columbia University Medical Center (CUMC) in the laboratory of Dr. Richard Sloan.

This document provides additional information about the protocol and variables included in the data file organized as follows: Overview of the protocol; description of the measures; detailed outline of the protocol and data processing; description of variables included in the data file, particularly the primary outcomes of interest and naming conventions.

Overview of Protocol

Purpose. The psychophysiology protocol in the MIDUS Biomarker project is a widely used laboratory based, experimental procedure designed to measure cardiovascular reactivity to and recovery from stress.

Procedure. During the protocol, participants' physiological outcomes are measured during a seated, resting baseline period followed by two cognitive/psychological stressor tasks, also in a seated position. The cognitive tasks are a mental arithmetic task called MATH (Morgan And Turner Hewitt arithmetic; Turner et al, 1986) and a Stroop color–word matching task. After each cognitive stress task, participants undergo another seated, resting period to assess physiological recovery to stress. The last period in the procedure is an orthostatic stressor. Participants move from a seated to standing position and remain standing for several minutes.

Physiological signal collection. Cardiovascular reactivity is assessed via continuous measurement of the electrocardiogram (ECG) and blood pressure measured at the finger and corrected to brachial artery standards. The beat-to-beat ECG and BP waveforms are then analyzed to calculate heart rate, several indices of heart rate variability (HRV), systolic and diastolic BP, and indices of BP variability (BPV). Heart rate variability is operationalized as variability in the series of intervals between consecutive R waves (the first upward deflection of the electrocardiogram following the Q-wave, arising from ventricular depolarization) (Figure 1). In addition, reactivity of the Hypothalamic Adrenal Pituitary (HPA) axis is measured via collection of saliva samples for cortisol assay.

Theory and method. Throughout this guide, relevant references are cited to provide investigators information about the methodology used in this protocol. We offer the following references, for an introduction and review of cardiac psychophysiology, and the types of questions that can be investigated with this type of protocol (Carney, Freedland, & Veith, 2005; Gorman & Sloan, 2000; Shcheslavskaya, Burg et al., (2010); Sloan, McCreath et al., 2007). A recent special issue of *Biological Psychology* (Allen & Chambers, 2007) on cardiac vagal control is a good resource. Investigators are encouraged to review the literature in more depth. Relevant key words for literature searches include: heart rate reactivity, heart rate variability (also referred to as “RR interval variability” or as a related measure, “respiratory sinus arrhythmia”), stress reactivity, and stress recovery.

The cognitive/psychosocial measures are described more fully in the next section. Details about the physiological measures and how they are processed appear later in this section after the comprehensive description of the protocol.

Psychological Stressor Measures

Stroop Color-Word Task. In this modified version of the Stroop task, one of four color name words (blue, green, yellow or red) is presented on a computer screen in a font color that is either congruent or incongruent with the name. The color name stimulus appears on screen, and participants press one of four keys on a keypad corresponding to the color of the letters in the word, not the color name. To roughly standardize the stressfulness of the task, the rate of presentation of the stimuli varies as a function of task performance. Greater accuracy leads to a more rapid presentation rate. Poorer accuracy leads to a slower rate. Overall, participants achieve an accuracy of about 67%. Response data for each trial, including the stimulus features, response latency, and response accuracy are stored in a file for later analysis.

Mental Arithmetic.

The Morgan And Turner Hewitt (MATH) task is a computer- administered mental arithmetic task designed for

use as a psychological stressor in laboratory studies of cardiovascular reactivity (Turner, Hewitt et al., 1986; Turner, Sims, Carroll, Morgan, & Hewitt, 1987). Task problems involve the addition or subtraction of two numbers. Problem difficulty can vary across five levels, ranging from problems of 1-digit \pm 1-digit numbers (level 1) to 3-digit \pm 3-digit numbers (level 5). The task always begins at level 3; difficulty level thereafter is determined at each trial by response accuracy on the previous trial. Correct responses were followed by one step up in difficulty, or if already at level 5, difficulty remains at level 5 until an incorrect response. Incorrect responses were followed by one step down in difficulty, or if already at level 1, difficulty remains at level 1 until a correct response.

Each trial consists of three elements presented on screen in succession. First, one math problem is presented for 2.0 sec. Then, the word 'Equals' appears alone on screen for 1.5 sec, giving the participant more processing time. A solution to the problem then appears for up to 1.0 sec, during which the participant presses one of two keys on a keypad to indicate whether the presented solution to the problem is correct or not. The next trial is presented as soon as a response key was pressed. Failure to respond within the one-second solution screen is recorded as an incorrect response, with a response time of 1.0 sec, and the next trial is presented.

Trials continued for the full duration of the mental arithmetic protocol period; total number of trials varies based on the participant's response times. The ratio of addition to subtraction problems is 3:7. The ratio of correct to incorrect problem solutions presented on screen is 1:1. Response data, including problem content, level, and response time and accuracy, are collected for each trial and stored in a file for later analysis. For MIDUS, the original task specifications by Turner et al. (Turner et al., 1986) were modified to extend the task length from 4 to 6 minutes.

Psychophysiology Protocol Description

The following is a detailed description of the data collection protocol, including equipment setup, protocol order, and data processing.

Protocol Flowsheet.

MIDUS staff who conducted this protocol used a data collection form called the psychophysiology flowsheet. A copy of this form with variable names inserted is in Section B (above).

The first two pages included questions about handedness, physical characteristics of the participant, and other factors that may influence experimental outcomes (e.g. consumption of caffeine, nicotine etc.), as well as a template of the protocol order. A more detailed version of this template appears in Table 1 below. The remainder of the flowsheet contains a more complete description of the protocol, instructions to research staff, instructions to participants, descriptions of the stress tasks, etc. Throughout the protocol, staff were instructed to record information at designated locations on flowsheet. This information as well as responses to the items at the beginning of the flowsheet were data-entered and included in the MIDUS Refresher Biomarker data file just prior to the physiology data described below and the variables begin with RA4V. Details about saliva sample collection are included in Section A above.

Monitoring Device Setup.

Electrocardiograph (ECG) electrodes were placed on the left and right shoulders, and in the left lower quadrant. Stretch bands were placed around the participant's chest and abdomen to measure respiration. A Finometer blood pressure cuff was placed on the middle finger of the non-dominant hand, and a Finometer blood pressure arm cuff was placed on the upper arm on the same side as the finger cuff. The participant was then seated and a numeric keypad, for responding to the stress tasks, was secured in a comfortable position relative to the dominant hand. The monitoring devices were then calibrated in the seated position.

Table 1. Detailed outline of psychophysiology protocol.

Protocol Order	Period	Period Description	Participant Activity	Period Duration (sec)	# Epochs
1		Saliva Sample #1 (upon arrival for session)	Salivette		
2		Saliva Sample #2 (Baseline)	Salivette		
3		Setup and Instructions	Practice stressor tasks	varies	
4	C1	Calibration: Finometer square wave & enter participant data.		≤ 600	
5	N1	Nada1: check physio signals' quality; reposition sensors if needed.		≤ 600	
6	C2	Calibration, seated: 1. Calibrate Finometer finger cuff to brachial artery (arm) BP. 2. Calibrate respiration volume.	1. Standard arm BP cuff inflates & deflates to auto-correct finger cuff BP readings. 2. 6 breaths into fixed volume 800ml plastic bag.	≤ 280	
7	B1	Baseline 1	Seated, quiet.	660	2
8	M1	Cognitive stressor 1	MATHTurner ^a	360 ^b	1
9	R1	Recovery 1		360	1
10	S1	Cognitive stressor 2	Stroop ^a	360	1
11	R2	Recovery 2		360	1
12	N2	Saliva Sample #3(After cognitive stressors & recovery) (Nada2 period)	Salivette	≤ 120	
13	N3	Nada3: Transition to standing position.		≤ 600	
14	C3	Calibration, standing: 1. Calibrate Finometer finger cuff to brachial artery (arm) BP. 2. Calibrate respiration volume.	1. Standard arm BP cuff inflates & deflates to auto-correct finger cuff BP readings. 2. 6 breaths into fixed volume 800ml plastic bag.	≤ 280	
15	U1	Physical (orthostatic) stressor: Standing upright		360	1
16		Saliva Sample #4 (After orthostatic stressor)	Salivette		
17		Saliva Sample #5 (30 min after Sample #4)	Salivette		

Notes:

^a The order of the Math and Stroop tasks was automatically and randomly selected at the time of data collection. Thus, each task was either in position 5 or 7 of the protocol for each session. A Task presentation order variable is included in the psychophysiology flowsheet data set.

Protocol Order.

The general protocol order was as follows (*details are in Table 1*): seated baseline (11 minutes); psychological stress task 1 (mental arithmetic or Stroop task - 6 minutes); recovery 1 (6 minutes); psychological stress task 2 (mental arithmetic or Stroop task - 6 minutes); recovery 2 (6 minutes); orthostatic stressor (standing/upright) (6 minutes). No recovery data were collected after exposure to the orthostatic stressor. Participants were instructed to remain silent throughout the procedures. After the second recovery period, participants were assisted in moving to a standing position. The monitoring devices were recalibrated, then the orthostatic stress period began.

Data Processing Criteria

The physiological monitoring equipment (ECG, Finometer, Inductotrace respirometer) ran continuously throughout the protocol and produced raw waveform data. These raw data were processed according to standardized algorithms (Task Force, 1996) to create variables (see [Key Variables](#)) that can be used in analyses. Analytic data are provided in MIDUS by **period** and by **epoch** within each period. The MIDUS Refresher Biomarker (P4) data includes one set of data from the psychophysiology session, which uses 300 second epochs of data. This section defines these terms and describes the criteria used to select raw physiological waveform data for processing to generate key variables.

Periods.

The protocol was divided into periods based on experimental conditions and participants' activity. Physiological outcome variables are computed separately for each protocol period and are identified by Period, as specified in Table 1, in the data sets.

Periods in **BLUE** font in Table 1 represent data included in the MIDUS data set for analyses.

Other periods represent interim periods used for calibrating equipment and other purposes not relevant to hypothesis testing. Raw physiological waveform data from the interim periods are preserved at the CUMC site but are not analyzed.

Each period name as shown in Table 1 is part of the variable names for all data from that period.

Data Epochs for Analysis

Within each protocol period, data were analyzed in specified epochs of time, based on different criteria and different types of research questions.

Data were analyzed with a specified 300 sec epoch duration.

For the 11 min baseline period, we attempt to provide 2 epochs of 300 sec each. Cases with unscorable intervals of data may have only 1 Baseline epoch (due to noisy signal). The stress tasks, recovery periods and upright stressor were all 6 min periods. For these, one epoch of data is included in the data set.

Salivary Cortisol Samples: Table 1 indicates the order and timing of the saliva samples. In the data set, saliva sample numbers (corresponding to numbers in Table 1) indicate the specific protocol context of the sample. For example, sample #4 was collected after the orthostatic stressor period. If sample #4 has a missing value in the data set, it means that participant does not have cortisol data related to orthostatic stress reactivity (period U1). Likely, it also indicates that the participant did not complete the U1 period in the protocol.

Physiological Measures

The following sections describe the 3 sets of physiological measures (HRV, BP, Respiration) and conventions

used in naming variables. Across all three sets of measures, the following characters are used to represent the indicated period when they were obtained during the protocol:

*Key to **Period** and **Epoch** indicators in variable names:*

B1	=	BASELINE period, epoch1 (of usually 2 epochs)
B2	=	BASELINE period, epoch2 (of usually 2 epochs)
M1	=	MATH period, epoch1 (of 1 epochs)
S1	=	STROOP period, epoch1 (of 1 epoch)
R1	=	first RECOVERY period (1 epoch)
R2	=	second RECOVERY period (1 epoch)
U1	=	UPRIGHT period (1 epoch)

A. Electrocardiogram Measures: Heart Rate, Heart Rate Variability

Acquisition and Processing of ECG Signals. Beat-to-beat analog ECG signals were collected then digitized at a sampling rate of 500 Hz by a 16-bit National Instruments analog-to-digital (A/D) board installed in a microcomputer. ECG waveforms were submitted to proprietary event detection software to identify R waves. Following established procedures, (Berntson, Quigley, Lang, & Boysen, 1990; Dykes, Ahmann et al., 1986), research staff visually reviewed all ECG waveforms to correct interactively any software errors in identifying normal R waves. The resulting series of normal RR intervals was used to calculate the cardiac variables heart rate (HR) and several standard indices of HRV.

Heart rate. Heart rate is calculated as an average of all valid RR intervals for a specified length of time. HR data in the MIDUS data set have been converted from RR interval units (milliseconds) to beats per minute units.

Heart rate variability. Time domain indices of RR interval variability include the standard deviation of RR intervals (SDRR) and the root mean squared successive differences (rMSSD). Frequency domain measures include spectral power in the low (0.04-0.15 Hz (LF-HRV)) and high (0.15-0.50 Hz (HF-HRV)) frequency bands. The spectra of RR interval series were calculated using an interval method for computing Fourier transforms similar to that described by DeBoer, et al. (DeBoer, Karemaker, & Strackee, 1984). Prior to computing Fourier transforms, the mean of the RR interval series is subtracted from each value in the series and the series then is filtered using a Hanning window (Harris, 1978) and the power, i.e., variance (in msec²), over the LF and HF bands is summed. Estimates of spectral power are adjusted to account for attenuation produced by this filter (Harris, 1978).

A. Key ECG Variables and Naming Conventions

a. Key Variables.

The key cardiac variables from the psychophysiology session used by CUMC investigators are listed below. These output variables are standardized based on conventions for measuring heart rate and heart rate variability parameters (Task Force, 1996) :

HR: Average heart rate, beats per minute units

SDRR: Standard deviation of RR intervals, msec units

rMSSD: Root mean squared successive differences, msec units

LF_HRV: Low frequency RR interval variability, bandwidth 0.04-0.15 Hz, msec² units

HF_HRV: High frequency RR interval variability, bandwidth 0.15-0.40 Hz, msec² units The data file includes both original and log transformed versions of all HRV variables (the last 4 variables

listed above) for each period, and each epoch within a given period, along with variables indicating Epoch Duration (secs) and the Number of R-R intervals analyzed in each epoch.

Note: the CUMC team always uses log-transformed versions of the variables (natural logarithm), a standard practice in HRV research, due to reliable skew in their distributions.

b. Variable Naming Conventions.

SPSS Variable Labels and Value Labels are included in the data set. Per MIDUS naming conventions, the variables derived from the ECG physiological signal have the following structure.

RA4Vaa bbb

Where:

RA4 = MIDUS Refresher (RA), Project 4 (4), per MIDUS conventions

V = MIDUS Code letter assigned to Psychophysiology Session data

The remaining characters are determined by the data. Since a common set of variables is generated for each period, the following conventions are applied to the key cardiac variables.

aa = Period Name and Epoch Number (e.g. B1, R2) corresponding to the conditions and tasks presented during psychophysiology session (see column 2 in Table 1). Most periods are represented by one character in the variable name. There were 2 Recovery periods, so they are shown as R1 or R2 to indicate their sequence in the protocol. The Baseline period was long enough to produce 2 epochs of physiological data. In these periods the second “a” character is used to indicate Epoch Number (1 or 2).

bbb = The final 2 or 3 characters identify the key outcome variables, described above, using the following abbreviations.

DU	Epoch duration in seconds
BEG	Begin/start time of epoch in elapsed seconds (from beginning of period)
END	Stop time of epoch in elapsed seconds (from beginning of period)
NU	Number of R-R intervals analyzed in epoch
HR	HR, Avg heart rate, beats per minute
SD	SDRR, standard deviation of R-R intervals, milliseconds
LSD	natural log of SDRR
RM	RMSSD, root mean square successive RR differences, milliseconds
LRM	natural log of RMSSD
LF	LF-HRV, low freq RR interval variability, bandwidth 0.04-0.15 Hz, msec ²
LLF	natural log of LF-HRV
HF	HF-HRV, high freq RR interval variability, bandwidth 0.15-0.40 Hz, msec ²
LHF	natural log of HF-HRV

Example Variable Names:

[RA4VB2HF] = Baseline period; Epoch 2; High Freq. HRV (not log transformed)

[RA4VR1LRM] = 1st Recovery period (R1); Epoch 1 (only 1 epoch for recovery periods); log-transformed rMSSD (Root mean squared successive differences)


B. Blood Pressure (BP) and Blood Pressure Variability (BPV): Overview

The psychophysiology protocol collected noninvasive measures of beat-to-beat blood pressure (BP) and blood

pressure variability (BPV). In humans, BPV measured in this way oscillates at the same frequencies as heart rate, i.e., between 0.003 – 0.50 Hz. Oscillations in the respiratory spectral range (0.15-0.50 Hz), termed “high frequency” (HF) are the mechanical product of respiration-induced intrathoracic pressure changes and, thus, are of little interest in psychophysiological research. BPV in the low frequency (LF) range (0.04-0.15 Hz), on the other hand, has been of considerable interest and may have physiologic significance, with evidence suggesting that it represents a central sympathetic oscillator, a resonance phenomenon, or the feedforward effects of heart rate variability (HRV) .

1. Scoring and Analysis of the BP Waveforms.

The beat-to-beat BP waveforms collected in this protocol go through an intensive process to identify physiologically valid data and use them to compute standard measures of BP and BPV. Systolic peaks and diastolic troughs are identified using custom-written software (see Figure below), resulting in 2 time series: the systolic and diastolic BP events. These time series are then visually reviewed by research staff trained to distinguish physiologically valid BP waveforms from artifact. Non-valid signals

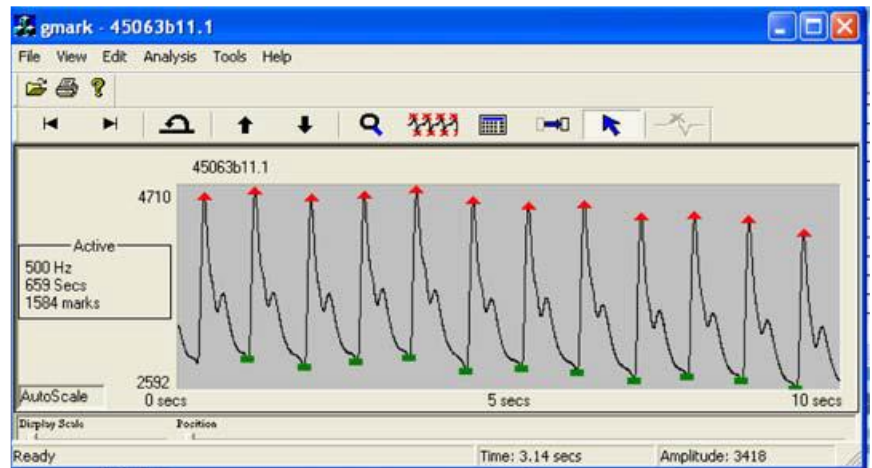

are corrected if possible via an interpolation algorithm then submitted to Fourier-based spectral analysis similar to that described by DeBoer et al. (1984), yielding estimates of power in the low frequency (LF, 0.04-0.15 Hz) band. Prior to computing Fourier transforms, the mean of the BP series is subtracted from each value in the series. The series is filtered using a Hanning window and the power over the LF band is summed. Estimates of spectral power are adjusted to account for attenuation produced by this filter (Harris, 1978).

2. Blood Pressure Variable Naming Conventions.

SPSS Variable Labels and Value Labels are included in the data set. Per MIDUS naming conventions the psychophysiology variable names have the following structure.

BP and BPV Root Variable Naming Scheme: **RA4Vbppcc**

Where:

RA4 = MIDUS Refresher (A), Project 4 (4),

V = MIDUS Code letter assigned to Psychophysiology Session data

The remaining characters are determined by the data type. Since a common set of variables is generated for all protocol periods, the following conventions are applied to all variables derived from the continuous BP waveforms measured during the protocol.

b = **BP component**; 1 character, identifies from which component of the blood pressure waveform the variable is derived.

S = SYSTOLIC BP (peak of BP waveform)

D = DIASTOLIC BP (trough of BP waveform)

pp = **Protocol Period and Data Epoch**; a one- or two-character sequence indicates the protocol period and sequential data epoch within the period represented by the variable. See the key to the naming scheme above .

cc = **Construct** being measured; 2 characters indicate the construct. These include the **key variables**, as well as **ancillary variables** that describe several aspects of the data collected during each period. These 2 sets of variables are described in the next section, along with measures of BaroReflex Sensitivity.

C. BP and BPV Variables

The data file contains three sets of variables in the following order: data quality codes, key and ancillary BP and BPV variables, and measures of BaroReflex Sensitivity (BRS). These sets of variables are described below beginning with the BP and BPV variables, then the data quality codes, and ending with BRS. Throughout the following, consistent with emerging nomenclature in experimental research, distinctions are made between key variables and ancillary variables. The key variables are the primary analytic variables derived from or computed using the physiological data. The ancillary variables include other descriptive or contextual data that users may want to consider when selecting variables for inclusion in analysis.

1. Key BP and BPV Variables.

- The key parameters derived from the blood pressure waveforms during each period, separately, of the psychophysiology protocol are listed below. These output variables are computed based on conventions for these constructs.
- Note that the data set includes **both a Systolic BP and Diastolic BP value for all the key and ancillary variables** listed below. That is, the systolic peak and diastolic trough of the BP waveforms are analyzed separately to produce 2 series of these key and ancillary variables.
- The data file includes both untransformed and natural log transformed versions of the LF-BPV variables because these variables reliably demonstrate positively skewed distributions. It is standard practice in the psychophysiology literature to log-transform these variables and use the transformed values for data analysis.

Var Name component	Construct	Measurement units
BP	Blood pressure	mmHg ^(a)
BS	Standard Deviation of blood pressure	mmHg ^(a)
LV	Low frequency blood pressure variability (LF-BPV), in spectral bandwidth from 0.04-0.15 Hz	msec ^{2(b)}
LL	Natural logarithm of LF-BPV value (transformation to correct for positively skewed distribution)	msec ^{2(b)}

Notes: ^(a) millimeters of mercury; ^(b) milliseconds squared

2. BP and BPV Ancillary Descriptive Variables.

- We provide these variables to give investigators descriptive information about the timing and duration of the physiological signals from which the computed key variables were derived.
- All data included in the data set can be considered valid for most purposes. Data considered nonvalid for use, based on the most common methods in the psychophysiology literature on blood pressure and BPV, have already been omitted from this data set.
- Nevertheless, some investigators, depending on the research questions, may choose to omit some periods of data based on information in these descriptive variables.

Var Name component	Construct	Measurement units
DU	Epoch duration of valid data used to compute key variables for this epoch	sec ^(a)
ST	start time of data epoch	elapsed sec
ET	end time of data epoch,	elapsed sec
MT	total recorded time during the period;	sec
IT	Total # BP-BP intervals recorded in this period	count
IU	# BP-BP intervals in this period that are valid and used to compute key variables	count
IA	average length of BP-BP interval	sec

Notes: ^(a) seconds

In the data file all the diastolic variables are listed first followed by all the systolic variables. Per MIDUS convention, ancillary variables are listed first followed by the key variables for each period.

EXAMPLE:

An example of the complete set of variable names for **one protocol period** is shown below. This example shows the list of all key and ancillary variables derived from *Diastolic BP data during the Math stressor task protocol period*.

RA4VDM1DU	Diastolic, Math Epoch1 duration, sec
RA4VDM1ST	Diastolic, Math Epoch1 start time, elapsed sec
RA4VDM1ET	Diastolic, Math Epoch1 end time, elapsed sec
RA4VDM1MT	Diastolic, Math Epoch1 total recorded time, sec
RA4VDM1IT	Diastolic, Math Epoch1 # B-B intervals total
RA4VDM1IU	Diastolic, Math Epoch1 # B-B intervals used
RA4VDM1IA	Diastolic, Math Epoch1 average length of B-B interval, sec
RA4VDM1BS	Diastolic, Math Epoch1, std dev of DBP, mmhg
RA4VDM1BP	Diastolic, Math Epoch1, DBP, mmhg
RA4VDM1LV	Diastolic, Math Epoch1, low freq BP variability 0.04-0.15 Hz, msec squared
RA4VDM1LL	Diastolic, Math Epoch1, natural log of LF BPV, msec squared

3. Baroreflex Sensitivity (BRS).

Blood pressure regulation – or dysregulation – is a crucial physiological function that can affect both acute and chronic health outcomes. Though blood pressure is regulated by several mechanisms, one of the most important is the baroreflex. The vascular system is equipped with stretch receptors called baroreceptors in the aortic arch and carotid arteries. BP changes on a beat-to-beat basis as blood is ejected from the heart, resulting in size fluctuations in arterial walls. The aortic and carotid baroreceptors signal these changes via the vagus nerve (cranial X) and carotid sinus nerve (to cranial IX) respectively. A homeostatic negative feedback reflex ensues, characterized by modulations of cardiac and vasomotor activity, which then result in changes in heart rate and vascular tone to offset changes in blood pressure (Hughson, et al., 1993; Voss, et al., 1999).

If the baroreflex is functioning properly, an increase in blood pressure results in a lengthening of the intervals between heartbeats i.e. R-wave to R-wave intervals (RRI). This RRI change allows BP to fall more between the next ventricular contractions than between shorter RRIs. Likewise, when BP decreases, baroreceptor feedback results in shorter RRIs.

In the MIDUS psychophysiology data, we can assess the cardiac arm of the baroreflex, as described above, because the aortic receptors signal via the vagus nerve, which also regulates heart rate variability.

This data set includes measures of baroreflex sensitivity, an index of how readily the heart responds to changes in BP. Reduced baroreflex sensitivity leads to increased transient changes in BP that may have acute or chronic ramifications.

a. BRS Indexes:

Traditionally, BRS has been measured as the RR interval (RRI) change corresponding to changes in BP produced by bolus injections of vasoactive drugs. Once regarded as the “gold standard,” most senior researchers in the field now believe that pharmacologic approaches are not free of problems and that noninvasive methods of BRS analysis are equally valid.

Accordingly, we use a sequence method approach. Series of **3 or more increasing or decreasing systolic BP values** constitute a “sequence” that likely represents engagement of the baroreflex. Such sequences are identified throughout each protocol period and used to compute BRS as described below.

b. Sequence Method Calculation steps: based on method described by Voss, et al. (1999), except that we used only the single sequence component of the method they described.


1. Identify corresponding BP and ECG data from each protocol period.
2. Identify and use only the Systolic BP values; remove Diastolic BP values
3. Compute RR time intervals from ECG data.
4. Remove Data occurring during Bad Intervals from both SBPs and RRs
5. Detect Sequences of 3 or more consecutive SBP increases or decreases
6. Match SBPs of identified Sequences with corresponding RR-Intervals
7. Use linear Regression to Estimate the Slope of Each Identified Sequence
8. Filter out unwanted Sequences (usually these are due to data scoring artifacts or errors)
 - a. Those with 2 BPs within 1 RR Interval,
 - b. Those with any RR Interval >2sec
9. Take mean slope (beta) of remaining Sequences

c. Baroreflex Sensitivity Root Variable Naming Scheme: **RA4Vppccc**

RA4V = MIDUS conventions, as stated earlier

pp = **Protocol Period**; two characters, using same characters as for the BP/BPV variable names. However, unlike for the BP/BPV variables, BRS values were **not** computed on epochs within protocol periods. Thus, there is a single set of BRS outcomes for each period. Investigators who want to compare periods with 2 epochs of BP/BPV outcomes with corresponding BRS outcomes can compute the mean of the 2 BP/BPV epochs.

ccc = **Construct** being measured; 3 characters indicate the construct. These include **key and ancillary variables** that describe several aspects of the data collected during each period. See tables below.

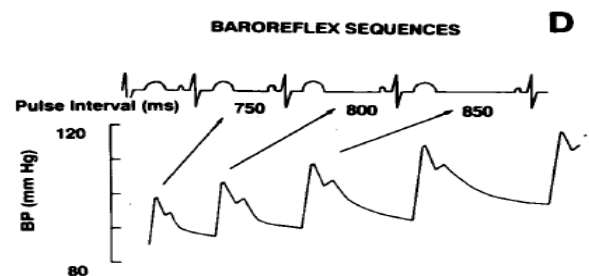

Key BRS Variables

As indicated by the last step in the sequence method calculation steps listed above, the primary index of BRS is computed as a **beta coefficient** i.e. the mean slope of the regression line produced by all sequences in each protocol period.

Var Name component	Construct	Measurement units
TMB	Total Mean Beta for All Sequences	ms/mmHg

Ancillary BRS Variables.

Var Name component	Construct	Measurement units
TNS	Total Number of Sequences in the period	Count

EXAMPLE: Baroreflex Sensitivity Variables (Note, the key and ancillary variables are listed separately above, but in the data file the Mean Beta for a given measure is listed first followed by the corresponding Total Number).

An example of the complete set of BRS variable names for **one protocol period**, the Math Stressor task is shown below.

RA4VM1TMB	Math, Mean Beta for All Sequences, ms/mmHg
RA4VM1TNS	Math, Total Number of Sequences

Respiration Rate

Respiration was monitored continuously during the protocol by inductive plethysmography using the Portable Inductotrace system (Bio-logic Systems Corp.®, Mundelein, Illinois). This device uses two stretch bands, one around the chest and abdomen, which measure volume excursions during the respiratory cycle. Respiratory volume was calibrated using an 800 ml plastic bag fitted with a mouthpiece tube (see psychophysiology protocol for information about respiration monitoring procedures). Analog signals from the two bands were sampled at 20 Hz and digitized, then summed into a single waveform for analysis. Using proprietary event detection software, respiratory excursions were identified automatically then corrected based on visual inspection to produce a breath-by-breath time series.

Proprietary analysis software computed respiration rate in 60 sec epochs from the breath-by-breath time series. Criteria for retaining respiration events as valid breaths in computing respiration rate were: a minimum volume ≥ 100 ml; and respiration rate of 6-30 breaths/minute. Missing values among the respiration variables in the 60 sec data set generally are due to movement artifacts or other sources of noisy physiological signals that prevented analysis of valid respiration signals during that minute of the protocol.

The respiration rate data, in 60 sec epochs which parallel the epochs provided for the HR and HRV variables, are provided so investigators can explore other questions or types of analyses.

Respiration Rate Variable Naming Scheme: RA4Vccppp

RA4V = MIDUS conventions, as stated earlier

cc = construct being measured (Respiration Rate) same characters (RR) for all variables,

pp = Protocol Period; two to three characters, using same characters indicating protocol period as indicated above, the final character indicates the Minute in the period.

Respiration Rate Adjustment

Heart rate variability in the high frequency range is influenced by respiration (Allen, Chambers, & Towers, 2007; Grossman, Karemaker, & Wieling, 1991; Grossman & Taylor, 2007; Grossman, Wilhelm, & Spoerle, 2004). In the psychophysiology literature, including much of the work by Dr. Richard Sloan, the lead investigator for this MIDUS stress psychophysiology protocol, HRV parameters are often adjusted for respiration rate prior to hypothesis testing analyses. That said, in Dr. Sloan's team, we often discover that adjusting for respiration makes no difference in the outcomes of analyses as compared with using unadjusted HRV values.

Because of this common practice in psychophysiology research, we provide these respiration data for investigators.

HRV measures adjusted for respiration are unstandardized residual scores estimating the variance in HRV that cannot be explained by the effect of respiratory rate. To compute respiration adjusted HF-HRV (corresponding in this data to the variables ending in "HF" e.g. RA4VB2HF):

1. compute the within-subject mean of all 60 sec epochs of respiration for each protocol period
2. conduct within-subject univariate regression analyses using respiratory rate as a predictor of HF-HRV (Sloan et al., 2001), on a period-by-period basis. The reason for computing each period separately is to avoid having within-subject effects mix with between subject effects in the error term if you were to include all protocol periods of data in the same regression analysis.
3. Save the residuals from these regression models as the new values of HF-HRV.
4. Finally, conduct a natural log transformation on these adjusted values of HRV. As you see, the data set includes log transformed values of all HRV parameters (e.g. RA4VB2LHF) because they are all routinely positively skewed. Thus, if you use respiration adjusted values of HF-HRV, you'll need to do the log transformation on these adjusted values.

SECTION D

DATA QUALITY FILTER VARIABLES

Psychophysiology Data Quality Filter Variables

Overview

The experimental psychophysiology protocol, conducted in the morning of the second day of the Biomarker visit, included assessments of beat-to-beat electrocardiogram (ECG), respiration, and beat-to-beat blood pressure. The data file contains standard outcome measures derived from each of these assessments.. Details about these measures, as well as the protocol, can be found in the preceding section (Section C)

If you have not read Section C please do so before continuing.

To facilitate analysis of the psychophysiology data, the data quality filter variables have been created for each set of assessments. This section describes these filter variables, which can be found in the data file immediately preceding the set of standard outcome measures that they are linked to.

Variables pertaining to the overall session are described first, followed by the period specific variables. Each set contains variables indicating the status of the session or period along with variables indicating the reason for that status designation. These variables can be used individually or in combination to select cases for inclusion in analyses.

Filter Variables: Period Specific

The beat-to-beat ECG waveform data for each period in the psychophysiology protocol were analyzed according to standard procedures, which are described elsewhere in this document. The quality (i.e. fidelity) of the ECG signal captured during data collection affects whether the data can be analyzed accurately. In addition, certain physiological anomalies (e.g. cardiac arrhythmias), represent non-sinus node neural input to the heart. R-R intervals that are not generated from the sinoatrial node are omitted from analysis. As a result of these various sources of error in the data, short or long sections of data in a “complete” session may show up as missing data in the final data set.

Data Quality Codes for the ECG and BP Data.

We provide with these data a set of variables describing the quality of the physiological data collected, whether nonvalid data occurred that resulted in missing values for the HR, HRV, BP and BPV outcome variables, and the nature of the nonvalid data. These data quality variables are meant to help investigators understand why missing data might exist even though the stress reactivity protocol was administered to a participant. Depending on one’s research questions, these quality codes may, thus, help make decisions on whether to include or omit certain cases from analyses. The data quality coding scheme was developed by Dr. Richard Sloan’s team at Columbia University Medical Center.

One important point for understanding these quality codes is that we reviewed and corrected the BP waveform data in conjunction with the concurrent ECG waveform data. The reasoning behind this procedure is that BP is primarily determined by cardiac ventricular contraction. If the timing of the ventricular contractions i.e. the ECG R waves, or the shape of the ECG waveforms represents a cardiac arrhythmia that is not valid for computing HRV, for our purposes in this study the corresponding BP waves also are not valid for use in computing BP or BPV. The valid computation of heart rate variability parameters depends on using ECG waves that, based on visual inspection, represent ***sinus rhythm*** of the R-R intervals i.e. contractions neurally generated via parasympathetic innervation of the sinoatrial (SA) node in the heart. For our interests in the MIDUS psychophysiological stress challenge protocol, we want to use corresponding blood pressure waves that are produced from normal sinus rhythm cardiac activity. Thus, anomalies in the ECG data are often used selectively to omit corresponding BP data.

Note that these ratings are made on the data quality of the overall BP waveforms for each protocol period. Unlike the computed BP and BPV variables, there is ***not*** a separate set of quality ratings for systolic and diastolic components of the waveforms.

Key Variables: Data Quality Codes for each Period:

Variable name component	Definition	Coded Values	Code definitions
WQ: Waveform quality	Physiological signal fidelity. If signal quality poor, valid signals cannot be identified for use in computing outcome parameters.	1=Clean signal	Period recorded clearly with no interference or noise.
		2=Noise	Period contains areas where signal is obscured by noise.
		3=Missing data.	Period recorded per physio protocol flowsheet but no valid data obtained; usually a technician error in which only the BP monitor's calibration waveform was recorded rather than participant data signal.
		4=Period not run.	Period not recorded due to early termination of protocol (determined by laboratory technician), or equipment malfunction; or entire session not run.
QR: Quality Reason	These codes help explain the reason for the Waveform Quality rating code.	1=Clean Signal	Waveform events are clearly visible; signal is consistently clean throughout file.
		2=Interference	Signal is obscured in some areas by electrical interference ("noise").
		3=PP movement	Signal morphology is contorted in some areas, likely due to participant movement.
		4=Loss of signal	Signal is lost in some areas, likely due to equipment malfunction or intentional termination of period (determined by laboratory technician).
		5=Missing data	Period recorded per physio protocol flowsheet but no valid data obtained; usually a technician error in which only the BP monitor's calibration waveform was recorded rather than participant data signal..
	PH: Physiology	6=Period not run	Period not recorded due to early termination of protocol (determined by laboratory technician), or equipment malfunction; or entire session not run.
		1=Normal	Physiology does not contain any arrhythmias identifiable by visual scan.
		2=PVCs	Premature ventricular contractions (PVCs) are one commonly encountered kind of non-sinus cardiac rhythm. Visible concurrently in ECG and BP signals. Can be validly corrected via interpolation, but if PVCs are too frequent (>20% of a file), effect on BP pulse timing and magnitude make the BP signal non-valid for analysis.
		3=PACs	Premature atrial contractions (PACs) are one commonly encountered kind of non-sinus cardiac rhythm. Visible concurrently in ECG and BP signals. Cannot be corrected , only omitted from period, and if too frequent may thus prevent there being enough usable continuous data to produce a long enough epoch from which to compute BP/BPV measures.
		4=Other Non-Sinus	Any arrhythmia not visually identifiable as PVC or PAC. A recurring interruption of pulse timing or magnitude making the BP signal non-valid for analysis.
		5=Missing data	Period recorded per physio protocol flowsheet but no valid data obtained; usually a technician error in which only the BP monitor's calibration waveform was recorded rather than participant data signal.
		6=Period not run	Period not recorded due to early termination of protocol (determined by laboratory technician), or equipment malfunction; or entire session not run.

BI: Bad Interval	Was there an unscorable interval of data omitted from analysis? If so, how much of the period was affected?	1=None	No unscorable intervals in the period.
		2=Partial	Some unscorable intervals have been identified in the period and omitted from analysis.
		3=Whole	The whole period has been identified as unscorable and omitted from analysis.
		4=Missing data	Period recorded per physio protocol flowsheet but no valid data obtained; usually a technician error in which only the BP monitor's calibration waveform was recorded rather than participant data signal.
		5=Period not run	Period not recorded due to early termination of protocol (determined by laboratory technician), or equipment malfunction; or entire session not run.
NO: Notes	Data scorer's notes.		Text data format. No numeric codes.

EXAMPLE: Quality Codes Variables

An example of the complete set of quality code variable names for ***one protocol period*** is shown below. This example shows the variables for data quality during the Math stressor task period.

RA4VM1BPWQ	Math, BP Waveform Quality
RA4VM1BPQR	Math, BP Quality Reason
RA4VM1BPPH	Math, BP Physiology
RA4VM1BPBI	Math, BP Bad Interval
RA4VM1NOT	Math, Scoring Notes

SECTION E

REFERENCES

References

- Allen, J. J. B., & Chambers, A. S. (Eds.). (2007). *Biological Psychology* (Vol. 74).
- Berntson, G. G., Quigley, K. S., Lang, J. F., & Boysen, S. T. (1990). An approach to artifact identification: Application to heart period data. *Psychophysiology*, 27, 586-598.
- Carney, R. M., Freedland, K. E., & Veith, R. C. (2005). Depression, the autonomic nervous system, and coronary heart disease. *Psychosomatic Medicine*, 67 Suppl 1, S29-33.
- DeBoer, R. W., Karemaker, J. M., & Strackee, J. (1984). Comparing spectra of a series of point events, particularly for heart rate variability spectra. *IEEE Transactions in Biomedical Engineering, BME-31*, 384-387.
- Dykes, F. D., Ahmann, P. A., Baldzer, K., Carrigan, T. A., Kitney, R. I., & Giddens, D. P. (1986). Breath amplitude modulation of heart rate variability in normal full term neonates. *Pediatric Research*, 20, 301-308.
- Gorman, J. M., & Sloan, R. P. (2000). Heart rate variability in depressive and anxiety disorders. *American Heart Journal*, 140(4 Suppl), 77-83.
- Grossman, P., & Taylor, E. W. (2007). Toward understanding respiratory sinus arrhythmia: Relations to cardiac vagal tone, evolution and biobehavioral functions. *Biological Psychology*, 74(2), 263-285.
- Grossman, P., Karemaker, J., & Wieling, W. (1991). Prediction of tonic parasympathetic cardiac control using respiratory sinus arrhythmia: The need for respiratory control. *Psychophysiology*, 28(2), 201-216.
- Grossman, P., Wilhelm, F. H., & Spoerle, M. (2004). Respiratory sinus arrhythmia, cardiac vagal control, and daily activity. *American Journal of Physiology - Heart and Circulatory Physiology*, 287(2 56-2), H728-H734.
- Harris, F. J. (1978). On the use of windows for harmonic analysis with the discrete Fourier transform. *Proceedings of the IEEE*, 66, 51-83.
- Hughson, R. L., Quintin, L., Annat, G., Yamamoto, Y., & Gharib, C. (1993). Spontaneous baroreflex by sequence and power spectral methods in humans. *Clinical Physiology*, 13(6), 663-676.
- Shcheslavskaya, O. V., Burg, M. M., McKinley, P. S., Schwartz, J. E., Gerin, W., Ryff, C. D., Weinstein, M., Seeman, T. E., & Sloan, R. P. (2010). Heart rate recovery after cognitive challenge is preserved with age. *Psychosomatic Medicine*, 72(2), 128-133. doi: 10.1097/PSY.0b013e3181c94ca0. PMC: 2950633.
- Sloan, R. P., Bagiella, E., Shapiro, P. A., Kuhl, J. P., Chernikhova, D., Berg, J., et al. (2001). Hostility, gender, and cardiac autonomic control. *Psychosomatic Medicine*, 63(3), 434- 440.
- Sloan, R., McCreath, H., Tracey, K., Sidney, S., Liu, K., & Seeman, T. (2007). RR Interval Variability Is Inversely Related to Inflammatory Markers: The CARDIA Study. *Molecular Medicine*, 13(3-4), 178-184.
- Task Force of the European Society of Cardiology and the North American Society of Pacing and Electrophysiology. (1996). Heart rate variability: Standards of measurement, physiological interpretation and clinical use. *Circulation*, 93(5), 1043-1065.
- Turner, J. R., Hewitt, J. K., Morgan, R. K., Sims, J., Carroll, D., & Kelly, K. A. (1986). Graded mental arithmetic as an active psychological challenge. *Int J Psychophysiol*, 3(4), 307-309.
- Turner, J. R., Sims, J., Carroll, D., Morgan, R. K., & Hewitt, J. K. (1987). A comparative evaluation of heart rate reactivity during MATH and a standard mental arithmetic task. *International Journal of Psychophysiology*, 5(4), 301-303.
- Voss, A., Baier, V., Malberg, H., Wessel, N., Hasart, A., Reinsperger, F., Schirdewan, A., & Osterziel, K. J. (1999, September 26-29). Application of the dual sequence method to estimate the baroreflex sensitivity in patients with dilated cardiomyopathy. Paper presented at the 26th Annual Meeting of Computers in Cardiology, Hannover, Germany. DOI: [10.1109/CIC.1999.825946](https://doi.org/10.1109/CIC.1999.825946)